DCM044-12  
Ed. 02/2024

## T3 ELISA

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## T3 ELISA

for routine analysis

Direct immunoenzymatic determination of Triiodothyronine (T3) in human serum or plasma.

IVD

LOT

See external label

2°C 8°C

Σ = 96 tests

REF DKO044

### INTENDED USE

Immunoenzymatic colorimetric method for quantitative determination of Triiodothyronine (T3) concentration in human serum and plasma from an adult population.

T3 ELISA kit is intended for laboratory use only.

### 1. CLINICAL SIGNIFICANCE

Triiodothyronine (T3) is a tyrosine-based hormones produced by the thyroid gland. An important component in the synthesis is iodine.

Thyroxine-binding globulin (TGB) is the major carrier protein for circulating thyroid hormone. Only the free fraction of T<sub>3</sub> is biologically active, a very small fraction of the circulating hormone is unbound - T3 0.3%.

The thyronines act on the body to increase the basal metabolic rate, affect protein synthesis and increase the body's sensitivity to catecholamines (such as adrenaline) by permissiveness. The thyroid hormones are essential for proper development and differentiation of all cells of the human body. These hormones also regulate protein, fat, and carbohydrate metabolism, affecting how human cells use energetic compounds. Numerous physiological and pathological stimuli influence thyroid hormone synthesis.

Both excess and deficiency of thyroxine can cause disorders. Thyrotoxicosis or hyperthyroidism is the clinical syndrome caused by an excess of circulating free thyroxine, free triiodothyronine, or both. It is a common disorder that affects approximately 2% of women and 0.2% of men. Hyperthyroidism is the case where there is a deficiency of thyroxine.


Release of T3 from the thyroid gland into the blood is ~5-8 micrograms/day. In addition, approximately 22 mg/day of T3 is produced by nonthyroidal 5-monodeiodination of T4. T3 has a much faster turnover than T4 (T ~1 day, compared to ~6 days for T4) and has greater biological potency than T4. Several drugs are known to affect the binding of Triiodothyronine to the thyroid hormone carrier proteins or its metabolism to T3 and complicate the interpretation of T3 results. Circulating autoantibodies to T3 and hormone-binding inhibitors may interfere.

Measurement of total serum T3 concentrations is a standard and well-validated test of thyroid gland function.

### 2. PRINCIPLE

The T3 (antigen) in the sample competes with the antigenic T3 conjugated with horseradish peroxidase (HRP) for binding to the limited number of antibodies anti T3 coated on the microplate (solid phase) (the enzyme conjugate should have no measurable binding to serum proteins especially TBG and albumin).

The interaction is illustrated by the following equation:

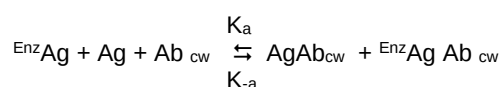Ab<sub>cw</sub>: monospecific immobilised antibody (constant quantity)

Ag: native antigen (variable quantity)

EnzAg: antigen conjugated to enzyme HRP (constant quantity)

Ag Ab<sub>cw</sub>: antigen-antibody complex

EnzAg Ab<sub>cw</sub>: antigen-HRP-antibody complex

K<sub>a</sub>: rate constant of association

K<sub>-a</sub>: rate constant of disassociation

K = k<sub>a</sub> / k<sub>-a</sub>: equilibrium constant

After incubation, the bound/free separation is performed by a simple solid-phase washing.

Then, the enzyme HRP in the bound-fraction reacts with the Substrate (H<sub>2</sub>O<sub>2</sub>) and the TMB Substrate and develops a blue color that changes into yellow when the Stop Solution (H<sub>2</sub>SO<sub>4</sub>) is added.

The colour intensity is inversely proportional to the T3 concentration of in the sample.

T3 concentration in the sample is calculated through a calibration curve.

### 3. REAGENTS, MATERIALS AND INSTRUMENTATION

#### 3.1. Reagents and materials supplied in the kit

##### 1. Calibrators (6 vials, 1 mL each)

|      |                   |
|------|-------------------|
| CAL0 | REF DCE002/4406-0 |
| CAL1 | REF DCE002/4407-0 |
| CAL2 | REF DCE002/4408-0 |
| CAL3 | REF DCE002/4409-0 |
| CAL4 | REF DCE002/4410-0 |
| CAL5 | REF DCE002/4411-0 |

##### 2. Conjugate Buffer (1 vial, 12.5 mL)

Buffer with preservative and binding proteins inhibitors

REF DCE002/4401-0

##### 3. Conjugate (1 vial, 1.4 mL)

T3 conjugated with horseradish peroxidase (HRP)

REF DCE002/4402-0

##### 4. Coated Microplate (1 breakable microplate)

Anti T3 antibody adsorbed on microplate

REF DCE002/4403-0

##### 5. TMB Substrate (1 vial) 15 mL

H<sub>2</sub>O<sub>2</sub>-TMB (0.26 g/L) (avoid any skin contact)

ProClin <0.0015%

REF DCE004-0

6. **Stop Solution** (1 vial, 15 mL)  
Sulphuric acid 0.15 mol/L (avoid any skin contact)  
**REF DCE005-0**
7. **50X Conc. Wash Solution** (1 vial, 20mL)  
NaCl 45 g/L, Tween-20 55 g/L.  
ProClin >0.0015% **REF DCE006-0**

### 3.2. Reagents necessary not supplied

Distilled water.

### 3.3. Auxiliary materials and instrumentation

Automatic dispenser.

Microplates reader (450 nm, 620-630 nm)

### Note

Store the reagents at 2-8°C in the dark.

Open the bag of reagent 5 (Coated Microplate) only when it is at room temperature and close it immediately after use; once opened, the microplate is stable until the expiry date of kit.

## 4. WARNINGS

- This kit is intended for in vitro use by professional persons only. Not for internal or external use in Humans or Animals.
- Use appropriate personal protective equipment while working with the reagents provided.
- Follow Good Laboratory Practice (GLP) for handling blood products.
- Some reagents (wash solution) contain small amounts of ProClin™ 300 (>0.0015%, <0.06%) as preservative. Avoid the contact with skin or mucosa.
- Classificazione secondo il regolamento (CE) n. 1272/2008 [CLP]  
Sensibilizzazione cutanea, categoria 1

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- The TMB Substrate contains an irritant, which may be harmful if inhaled, ingested or absorbed through the skin. To prevent injury, avoid inhalation, ingestion or contact with skin and eyes.
- The Stop Solution consists of a diluted sulphuric acid solution. Sulphuric acid is poisonous and corrosive and can be toxic if ingested. To prevent chemical burns, avoid contact with skin and eyes.
- Avoid the exposure of reagent TMB/H<sub>2</sub>O<sub>2</sub> to directed sunlight, metals or oxidants. Do not freeze the solution.
- Not intended for newborn screening

## 5. PRECAUTIONS

- Please adhere strictly to the sequence of pipetting steps provided in this protocol. The performance data represented here were obtained using specific reagents listed in this Instruction For Use.
- All reagents should be stored refrigerated at 2-8°C in their original container. Any exceptions are clearly indicated. The reagents are stable until the expiry date when stored and handled as indicated.
- Allow all kit components and specimens to reach room temperature (22-28°C) and mix well prior to use.
- Do not interchange kit components from different lots. The expiry date printed on box and vials labels must be observed. Do not use any kit component beyond their expiry date.
- If you use automated equipment, the user has the responsibility to make sure that the kit has been appropriately validated for its intended use/purpose.
- The incomplete or inaccurate liquid removal from the wells could influence the assay precision and/or increase the background. To improve the performance of the kit on automatic systems, it is recommended to increase the number of washes.
- It is important that the time of reaction in each well is held constant for reproducible results. Pipetting of samples should not extend beyond ten minutes to avoid assay drift. If more than 10 minutes are needed, follow the same order of dispensation. If more than one plate is used, it is recommended to repeat the dose response curve in each plate
- Addition of the TMB Substrate solution initiates a kinetic reaction, which is terminated by the addition of the Stop Solution. Therefore, the TMB Substrate and the Stop Solution should be added in the same sequence to eliminate any time deviation during the reaction.
- Observe the guidelines for performing quality control in medical laboratories by assaying controls and/or pooled sera.
- Maximum precision is required for reconstitution and dispensation of the reagents.
- Samples microbiologically contaminated, highly lipemic, icteric or haemolysed should not be used in the assay.
- Plate readers measure vertically. Do not touch the bottom of the wells.
- Fresh disposable tips must be used when pipetting assay reagents including samples, calibrators and controls to mitigate the risk of carryover contamination. Failure to do so may lead to invalid results.

## 6. PROCEDURE

### 6.1. Preparation of the Calibrators (C<sub>0</sub>...C<sub>5</sub>)

The Calibrators are ready to use, are calibrated against Human Serum Reference for triiodothyronine and have approximate concentrations of:

|       | C <sub>0</sub> | C <sub>1</sub> | C <sub>2</sub> | C <sub>3</sub> | C <sub>4</sub> | C <sub>5</sub> |
|-------|----------------|----------------|----------------|----------------|----------------|----------------|
| ng/dL | 0              | 50             | 100            | 250            | 500            | 750            |

**Exact levels are given on the labels on a lot specific basis.** Once opened, the Calibrators are stable 6 months at 2-8°C.

### 6.2. Preparation of Diluted Conjugate

Dilute the T3 conjugate 1:11 with Total T3 Conjugate Buffer (i.e. dilute 160 µL of conjugate with 1.6 mL of conjugate buffer).

Stable 24 h at 2-8°C.

### 6.3. Preparation of Wash Solution

Dilute the contents of each vial of the "50X Conc. Wash Solution" with distilled water to a final volume of 1000 mL prior to use. For smaller volumes respect the 1:50 dilution ratio. The diluted wash solution is stable for 30 days at 2÷8°C.

### 6.4. Preparation of the Sample

The determination of T3 should be performed in human serum or plasma.

Specimens may be refrigerated at 2÷8°C (for a maximum period of 48 hours). If the specimens cannot be assayed within 48 hours, the samples may be stored at temperatures of -20°C for up to 30 days. Avoid repetitive freezing and thawing of samples. When assayed in duplicate, 0,10 mL of the specimen are required.

### 6.5. Procedure

- **Allow all reagents to reach room temperature (22-28°C) for at least 30 minutes.** At the end of the assay, store immediately the reagents at 2-8°C: avoid long exposure to room temperature.
- Unused coated microwell strips should be released securely in the foil pouch containing desiccant and stored at 2-8°C.
- To avoid potential microbial and/or chemical contamination, unused reagents should never be transferred into the original vials.
- As it is necessary to perform the determination in duplicate in order to improve accuracy of the test results, prepare two wells for each point of the calibration curve (C<sub>0</sub>-C<sub>5</sub>), two for each sample, one for Blank.

| Reagent                                                                                                                                                                                                                                                                                                                                                                                                                   | Calibrator | Sample/<br>Control | Blank  |
|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|------------|--------------------|--------|
| Sample                                                                                                                                                                                                                                                                                                                                                                                                                    |            | 50 µL              |        |
| Calibrator<br>C <sub>0</sub> -C <sub>5</sub>                                                                                                                                                                                                                                                                                                                                                                              | 50 µL      |                    |        |
| Diluted<br>Conjugate                                                                                                                                                                                                                                                                                                                                                                                                      | 100 µL     | 100 µL             |        |
| Incubate 1h at room temperature (22÷28°C).<br>Remove the content from each well, wash the wells 3 times with 300 µL of diluted Wash Solution.<br><b>Important note:</b> during each washing step, gently shake the plate for 5 seconds and remove excess solution by tapping the inverted plate on an absorbent paper towel.<br><b>Automatic washer:</b> if you use automated equipment, wash the wells at least 5 times. |            |                    |        |
| TMB<br>Substrate                                                                                                                                                                                                                                                                                                                                                                                                          | 100 µL     | 100 µL             | 100 µL |
| Incubate 15 minutes in the dark at room temperature (22÷28°C).                                                                                                                                                                                                                                                                                                                                                            |            |                    |        |
| Stop<br>Solution                                                                                                                                                                                                                                                                                                                                                                                                          | 100 µL     | 100 µL             | 100 µL |
| Shake the microplate gently.<br>Read the absorbance (E) at 450 nm against a reference wavelength of 620-630 nm or against Blank within 5 minutes.                                                                                                                                                                                                                                                                         |            |                    |        |

## 7. QUALITY CONTROL

Each laboratory should assay controls at levels in the hypothyroid, euthyroid and hyperthyroid range for monitoring assay performance. These controls should be treated as unknowns and values determined in every test procedure

performed. Quality control charts should be maintained to follow the performance of the supplied reagents. Pertinent statistical methods should be employed to ascertain trends. The individual laboratory should set acceptable assay performance limits. Other parameters that should be monitored include the 80, 50 and 20% intercepts of the Calibration curve for run-to-run reproducibility. In addition, maximum absorbance should be consistent with past experience. Significant deviation from established performance can indicate unnoticed change in experimental conditions or degradation of kit reagents. Fresh reagents should be used to determine the reason for the variations.

In order for the assay results to be considered valid the kit calibrators and control must fall within the specifications detailed in the lot specific certificate of analysis.

If a control is out of its specified range, the associated test results are invalid and samples must be retested.

## 8. RESULTS

### 8.1. Mean Absorbance

Calculate the mean of the absorbance (E<sub>m</sub>) for each point of the Calibration curve (C<sub>0</sub>-C<sub>5</sub>) and of each sample.

### 8.2. Calibration curve

Plot the mean value of absorbance (E<sub>m</sub>) of the Calibrators (C<sub>0</sub>-C<sub>5</sub>) against concentration. Draw the best-fit curve through the plotted points. (es: Four Parameter Logistic).

### 8.3. Calculation of Results

Interpolate the values of the samples on the Calibration curve to obtain the corresponding values of the concentrations expressed in ng/dL.

The conversion factor for the unit, sometimes reported in other systems is:

100 ng/dL = 1 ng/mL

## 9. REFERENCE VALUES

A study of euthyroid adult population was undertaken to determine expected values for the T3 ELISA kit.

|       | Mean (ng/dL) | SD   | Range (ng/dL) |
|-------|--------------|------|---------------|
| Value | 118,5        | 33,4 | 52-185        |

Please pay attention to the fact that the determination of a range of expected values for a "normal" population in a given method is dependent on many factors, such as specificity and sensitivity of the method used and type of population under investigation. Therefore each laboratory should consider the range given by the Manufacturer as a general indication and produce their own range of expected values based on the indigenous population where the laboratory works.

## 10. PERFORMANCE AND CHARACTERISTICS

### 10.1. Precision

#### 10.1.1. Intra Assay Variation

Within run variation was determined by replicate determination (20x) of two different control sera in one assay. The within assay variability is ≤ 10.7%.

#### 10.1.2. Inter Assay Variation

Between run variations was determined by replicate measurements of three different control sera in different lots of kit. The between assay variability is  $\leq 9.1\%$ .

#### 10.2. Accuracy

The recovery of 50 – 100 – 200 – 400 ng/dL of T3 added to sample gave an average value ( $\pm$ SD) of  $97.5\% \pm 4\%$  with reference to the original concentrations.

#### 10.3. Sensitivity

The lowest detectable concentration of T3 that can be distinguished from the Calibrator 0 is 5 ng/dL at the 95 % confidence limit.

#### 10.4. Specificity

The cross-reactivity of the triiodothyronine antibody to selected substances was evaluated by adding the interfering substance to a serum matrix at various concentrations. The cross-reactivity was calculated by deriving a ratio between doses of interfering substance to dose of Triiodothyronine needed to displace the same amount of tracer.

| Substance                 | Cross Reactivity | Conc.     |
|---------------------------|------------------|-----------|
| l-Triiodo-thyronine       | 1.0000           | -         |
| l-Thyroxine               | 0.01             | 10 µg/dL  |
| d-Thyroxine               | 0.0025           | 10 µg/dL  |
| d-Triiodo-thyronine       | 0.015            | 100 µg/dL |
| Monoiodo-Tyrosine         | N/D              | 100 µg/mL |
| Diiodo-Tyrosine           | N/D              | 100 µg/mL |
| Triiodothyroacetic Acid   | N/D              | 100 µg/mL |
| Tetraiodothyroacetic Acid | N/D              | 100 µg/mL |

#### 10.5. Correlation with RIA

Dia.Metra T3 ELISA was compared to another commercially available T3 assay. 140 serum samples were analysed according in both test systems.

The linear regression curve was calculated:

$$(T3 \text{ Dia.Metra}) = 0.990 \cdot (T3 \text{ RIA}) + 0.6025$$

$$r^2 = 0.982$$

### 11. WASTE MANAGEMENT

Reagents must be disposed off in accordance with local regulations.

All materials that have come into contact with samples and reagents must be disposed of in accordance with country, state and local regulations.

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